

HEK293/Human GPC3 ΔHS Stable Cell Line

Catalog No.	Size
CHEK-ATP212	$2 \times (1 \text{ vial contains } \sim 5 \times 10^{6} \text{ cells})$

• Description

The HEK293/Human GPC3 Δ HS Stable Cell Line was engineered to express the human GPC3 Δ HS lacking the heparan sulfate (HS) chains, which produced by introducing the point mutations (S495A and S509A) on the wild-type human GPC3 (Gene ID: 2719). Surface expression of human GPC3 Δ HS was confirmed by flow cytometry.

• Application

Useful for cell-based GPC3 ΔHS binding assay

• Cell Line Profile

Cell line	HEK293/Human GPC3 ΔHS Stable Cell Line	
Host Cell	HEK293	
Property	Adherent	
Complete Growth Medium	DMEM + 10% FBS	
Selection Marker	Puromycin (2 μg/mL)	
Incubation	37°C with 5% CO ₂	
Doubling Time	22-24 hours	
Transduction Technique	Lentivirus	



• Materials Required for Cell Culture

- DMEM medium (Gibco, Cat. No. 11965-092)
- Fetal bovine serum (CellMax, Cat. No. SA211.02)
- Puromycin (InvivoGen, Cat. No. ant-pr-5b)
- 0.25% Trypsin-EDTA (1X), Phenol Red (Gibco, Cat. No. 25200-056)
- Penicillin-Streptomycin (Gibco, Cat. No. 15140-122)
- Phosphate Buffered Saline (1X) (HyClone, Cat. No. SH30256.01)
- Complete Growth Medium: DMEM + 10% FBS, 1% P/S
- Culture Medium: DMEM + 10% FBS, Puromycin (2 μg/mL), 1% P/S
- Freeze Medium: 90% FBS, 10% (V/V) DMSO
- T-75 Culture flask (Corning, 430641)
- Cryogenic storage vials (SARSTEDT, 72.379.007)
- Thermostat water bath
- Centrifuge
- Luna cell counter (Logos Biosystems, LUNA-II)
- CO₂ Incubator (Thermo, 3111)
- Biological Safety Cabinet (Thermo, 1389)



• Recovery

- 1. Thaw the vial by gentle agitation in a 37°C water bath. To reduce the possibility of contamination, keep the cap out of the water. Thawing should be rapid (approximately 2 minutes).
- 2. Remove the vial from the water bath as soon as the contents are thawed, and decontaminate by spraying with 70% ethanol. All the operations from this point on should be carried out under strict aseptic conditions.
- 3. Transfer the vial contents to a centrifuge tube containing 4.0 mL complete growth medium and spin at approximately 1000 rpm for 5 minutes.
- 4. Resuspend cell pellet with 5 mL complete growth medium and transfer the cell suspension into T-75 flask containing 10-15 mL of pre-warmed complete growth medium.
- 5. Incubate at 37°C with 5% CO₂ incubator until the cells are ready to be split.

• Subculture

- 1. Remove and discard culture medium.
- 2. Wash the cells once with sterile PBS.
- 3. Add 2 mL of 0.25% trypsin to cell culture flask. Place the flask at 37°C for 2-3 minutes, until 90% of the cells have detached.
- 4. Add 6.0 to 8.0 mL of culture medium and aspirate cells by gently pipetting.
- 5. Add appropriate aliquots of the cell suspension to new culture vessel.
- 6. Incubate at 37℃ with 5% CO₂ incubator.

Subcultivation Ratio: A subcultivation ratio of 1:6 to 1:10 is recommended.

Medium Renewal: Every 2 to 3 days.

Note: After recovery for 1-2 generations with the complete growth medium not containing the selection marker, if the cell state is well, changing to the culture medium containing the selection marker.



• Cryopreservation

- 1. Remove and discard spent medium.
- 2. Detach cells from the cell culture flasks with 0.25% trypsin.
- 3. Centrifuge at 1000 rpm for 5 min at RT to pellet cells.
- 4. Resuspend the cell pellets with complete growth medium and count viable cells.
- 5. Centrifuge at 1000 rpm for 5 min at RT and resuspend cells in freezing medium to a concentration of 5×10^6 to 1×10^7 cells/mL.
- 6. Aliquot into cryogenic storage vials. Place vials in a programmable cooler or an insulated box placed in a –80°C freezer overnight, then transferring to liquid nitrogen storage.


• Storage

- Product format: Frozen
- Storage conditions: Liquid nitrogen immediately upon receipt



• Receptor Assay

using PE-labeled anti-human GPC3 antibody.

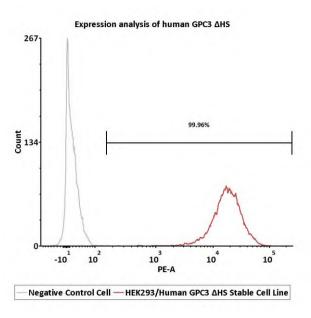

Catalog No.	Stable Cell Line	MFI for GPC3 ΔHS (PE)
NA	Negative Control Cell	10.65
CHEK-ATP212	HEK293/Human GPC3 AHS Stable Cell Line	17297.73

Fig1. Expression analysis of human GPC3 ΔHS on HEK293/Human GPC3 ΔHS Stable Cell Line by FACS. Cell surface staining was performed on HEK293/Human GPC3 ΔHS Stable Cell Line or negative control cell



• License Disclosure

This reporter cell is provided for research use only. This license does not permit you to share, distribute, sell, sublicense, or otherwise make this reporter cell available for use to other laboratories, departments, research institutions, hospitals, universities, or biotech companies. The license does not permit modification of this reporter cell in any way. Inappropriate use or distribution of this reporter cell will result in revocation of the license. Modifications of this cell line, transfer to another facility, or commercial use of the cells may require a separate license and additional fees. AcroBiosystems does not warrant the suitability of this reporter cell for any particular use, and does not accept any liability in connection with the handling or use of this reporter cell.

• Related Products

<u>Products</u>	Cat.No.	
CHO/Human Mesothelin Stable Cell Line Development Service	SCCHO-ATP120	
CHO/Human Glypican-3 (GPC3) Stable Line Development Service	SCCHO-ATP112	
CHO/Human STEAP1 Stable Cell Line Development Service	SCCHO-ATP121	
CHO/Human uPAR Stable Cell Line Development Service	SCCHO-ATP152	
CHO/Human c-MET Stable Cell Line Development Service	SCCHO-ATP141	
HEK293/Human ROR1 Stable Cell Line	CHEK-ATP084	
HEK293/Human Mesothelin Stable Cell Line	CHEK-ATP119	
HEK293/Human Glypican-3 (GPC3) Stable Cell Line	CHEK-ATP092	
HEK293/Human FOLR1 Stable Cell Line	CHEK-ATP091	
HEK293/Human DLL3 Stable Cell Line	CHEK-ATP090	
HEK293/Human TL1A Stable Cell Line	CHEK-ATP142	
HEK293/Human NAPI-IIb Stable Cell Line	CHEK-ATP116	
HEK293/Human Cadherin-6 Stable Cell Line	CHEK-ATP127	
HEK293/Human ENPP3 Stable Cell Line	CHEK-ATP122	
HEK293/Human B7-H4 Stable Cell Line	CHEK-ATP126	
HEK293/Human Cadherin-17 Stable Cell Line	CHEK-ATP173	
HEK293/Human TL1A Stable Cell Line	CHEK-ATP142	
MDCK/Mouse FCGRT-P2A-mGFP&B2M Cell Line Development Service	SCMDC-ATP196	