

Rabies virus Nucleoprotein ELISA Kit

Pack Size: 96 tests

Catalog Number: RAS-A195

IMPORTANT: Please carefully read this manual before performing your experiment.

For Research Use Only. Not For Use in Diagnostic or Therapeutic Procedures





INTENDED USE

This kit is developed for quantitative detection of Rabies virus Nucleoprotein in samples.

It is intended for research use only (RUO).

PRINCIPLE OF THE ASSAY

Rabies virus (RABV), scientific name Rabies lyssavirus, is a deadly neurotropic virus that causes rabies in humans and animals. Rabies virus has an extremely wide host range and its transmission most often occurs through the saliva of animals. Without intervention prior to disease progression, rabies has the highest case fatality of any infectious disease. RABV contains a single-stranded negative-sense R genome that encodes five structural proteins: nucleoprotein (N), phosphoprotein (P), matrix protein (M), glycoprotein (G), and R-dependent R polymerase (L). Among these viral proteins, the RABV glycoprotein (RABV-G) is a pivotal player mediating virus entry and the major target of neutralizing antibodies, thus a key factor for vaccine and drug design. A rapid and effective assay kit detecting the levels of RABV-N is urgently needed to accelerate the development of RABV vaccines.

This assay kit is used to measure the levels of Rabies virus Nucleoprotein by employing a standard sandwich-ELISA format. The microplate in the kit has been pre-coated with Anti-Nucleoprotein (RABV) Antibody. First add the standard samples provided in the kit and your samples to the plate, incubate and wash the wells. Then add the Biotin-Anti-Nucleoprotein (RABV) Antibody to the plate, incubate and wash the wells. Next add Streptavidin-HRP to the plate, incubate and wash the wells. Lastly load the substrate into the wells and monitor color development in proportion with the amount of Rabies virus Nucleoprotein protein present. The reaction is stopped by the addition of a stop solution and the intensity of the absorbance can be measured at 450 nm and 630 nm. The OD Value reflects the amount of Rabies virus Nucleoprotein bound.

MATERIALS PROVIDED

TABLE 1. MATERIALS PROVIDED

Catalog	Components	Size	Format	Storage		
Cuturog	Components	(96 tests)	1 01	Unopened	Opened	
RAS195-C01	Pre-coated Anti-Nucleoprotein (RABV) Antibody Microplate	1 plate	Solid	2-8°C	2-8°C	

2 / 10

US and Canada:
Asia and Pacific:

Tel: +1 800-810-0816

Tel: +86 400-682-2521

E-mail: order@acrobiosystems.com



RA195-EN.01

RAS195-C02	Nucleoprotein (RABV) Standard	15 μg	Powder	2-8°C	-70°C
RAS195-C03	Biotin-Anti- Nucleoprotein (RABV) Antibody	15 μg	Powder	2-8°C	-70°C
RAS195-C04	Streptavidin-HRP	15 μg	Powder	2-8°C, avoid light	-70°C, avoid light
RAS195-C05	10×Washing Buffer	50 mL	Liquid	2-8°C	2-8°C
RAS195-C06	2×Dilution Buffer	50 mL	Liquid	2-8°C	2-8°C
RAS195-C07	Substrate Solution	12 mL	Liquid	2-8°C, avoid light	2-8°C, avoid light
RAS195-C08	Stop Solution	7 mL	Liquid	2-8°C	2-8°C

REAGENTS/EQUIPMENT NEEDED BUT NOT SUPPLIED

Single or dual wavelength microplate reader with 450 nm and 630 nm filter;

Centrifuge;

37°C Incubator;

10 μL, 200 μL and 1000 μL precision pipettes;

 $10 \mu L$, $200 \mu L$ and $1000 \mu L$ pipette tips;

Multichannel pipettes;

Tubes;

Graduated cylinder to prepare Wash Solution;

Deionized or distilled water to dilute 10×Washing Buffer;

STORAGE

- 1. Unopened kit should be stored at 2°C-8°C upon receiving.
- 2. Find the expiration date on the outside packaging and do not use reagents past their expiration date.
- 3. The opened kit should be stored per components table. The shelf life is 30 days from the date of opening.

REAGENT PREPARATION

1. Bring all reagents and samples to room temperature (20°C-25°C) before use. If crystals have formed in buffer solution, place the sample in a 37 °C incubator until the crystals have completely dissolved and bring the solution back

3 / 10

US and Canada: Tel: +1 800-810-0816 Web: http://www.acrobiosystems.com
Asia and Pacific: Tel: +86 400-682-2521 E-mail: order@acrobiosystems.com



to room temperature before use.

2. Reconstitute the provided lyophilized materials to stock solutions with distilled, sterile water as recommended in Table 2 and place the materials for 15 to 30 minutes at room temperature with occasional gentle mixing. Avoid vigorous shaking. The reconstituted stock solutions should be stored at -70°C. It is recommended not to freeze-thaw more than 1 time, the packing specification shall not be less than 5 μg.

TABLE 2. RECONSTITUTION METHODS FOR 96 TESTS

ID	Components	Size	Stock Solution Con.	Reconstitution Buffer and Vol.
RAS195-C02	Nucleoprotein (RABV) Standard	15 μg	100 μg/mL	150 μL water
RAS195-C03	Biotin-Anti- Nucleoprotein (RABV) Antibody	15 μg	100 μg/mL	150 μL water
RAS195-C04	Streptavidin-HRP	15 μg	100 μg/mL	150 μL water

RECOMMENDED SAMPLE PREPARATION

1. Working Fluid Preparation

1.1 Preparation of 1×Washing Buffer:

Dilute 50 mL 10×Washing Buffer with ultrapure water/deionized water to 500 mL.

1.2 Preparation of 1×Dilution Buffer:

Dilute 50 mL 2×Dilution Buffer with 1×Washing Buffer to 100 mL.

1.3 Preparation of Biotin-Anti- Nucleoprotein (RABV) Antibody working fluid:

Dilute Biotin-Anti- Nucleoprotein (RABV) Antibody to $0.2 \mu g/mL$ with Dilution Buffer. Please prepare it for one-time use only.

1.4 Preparation of Streptavidin-HRP working fluid:

Dilute Streptavidin-HRP to $0.1 \mu g/mL$ with $1 \times Dilution$ Buffer. The prepared working fluid should avoid light. Please prepare it for one-time use only.

2. Preparation of Standard Curve

Make serial dilutions of the Nucleoprotein (RABV) as a Standard curve with Dilution Buffer as recommended in Figure 1.


4 / 10

US and Canada: Tel: +1 800-810-0816 Web: http://www.acrobiosystems.com
Asia and Pacific: Tel: +86 400-682-2521 E-mail: order@acrobiosystems.com





FIGURE 1. PREPARATION OF 1:1 SERIAL DILUTIONS OF THE Nucleoprotein (RABV)

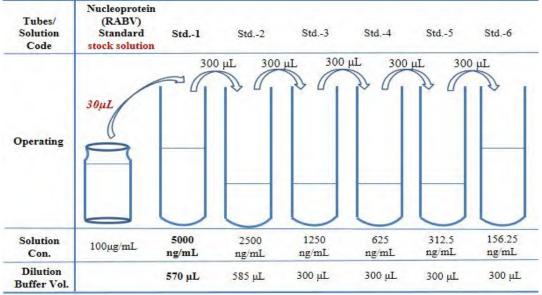

3. Add Samples

Add 100μL serially diluted **Nucleoprotein (RABV)** Standard curve and samples to each well. For blank Control wells, please add 100μL 1×Dilution Buffer. Seal the plate with microplate sealing film and incubate at room temperature for 1.0 hour.

Note: It is recommended to set double holes for samples and standard curves to be tested.

4. Washing

Remove the remaining solution by aspiration, add 300 µL of 1×Washing Buffer to each well, soak for 30s, remove any remaining 1×Washing Buffer: by aspirating or decanting, invert the plate and blot it against paper towels. Repeat the wash step above for three times.

5. Add Biotin-Anti-Nucleoprotein (RABV) Antibody

For all wells, add 100 µL Biotin-Anti-Nucleoprotein (RABV) Antibody (dilute to 0.2 µg/mL) working solution. Seal the plate with microplate sealing film and incubate at room temperature for 1.0 hour.

6. Washing

Repeat step 4.

7. Add Streptavidin-HRP

5 / 10

US and Canada: **Tel**: +1 800-810-0816

E-mail: order@acrobiosystems.com



ACTO*

For all wells, add 100 µL **Streptavidin-HRP** (**dilute to 0.1 µg/mL**) working solution. Please prepare it for one-time use only, avoid light. Seal the plate with microplate sealing film and incubate at room temperature for 1.0 hour.

8. Washing

Repeat step 4.

9. Substrate Reaction

Add 100 µL **Substrate Solution** to each well. Seal the plate with microplate sealing film and incubate at room temperature for 20 min, avoid light.

10. Termination

Add 50 µL Stop Solution to each well and tap the plate gently to allow thorough mixing.

Note: the color in the wells should change from blue to yellow.

11. Data Recording

Read the absorbance at 450 nm and 630 nm using UV/Vis microplate spectrophotometer within 5 minutes.

Note: To reduce the background noise, subtract the value read at $OD_{450 \text{ nm}}$ with the value read at $OD_{630 \text{ nm}}$.

CALCULATION OF RESULTS

- 1. Normal range of Standard curve: R²≥0.9900, detection range: 156.25-5000 ng/mL.
- 2. If the OD value of the sample to be tested is higher than the highest standard, the sample shall be diluted with dilution buffer and assay repeated.
- 3. To calibrate absorbance value obtained by the standard curve, the OD value of the sample to be measured is subtracted from the OD value of the blank control. The standard curve is plotted with the standard concentration as x-axis and the calibrated absorbance value as y-axis. Four parameters logistic are used to draw the standard curve and calculate the sample concentration.

PRECAUTIONS

- 1. This kit is for research use only and is not for use in diagnostic or therapeutic procedures.
- 2. The kit should be used according to the instructions.
- 3. Do not mix reagents from different lots.
- 4. Bring all reagents and samples to room temperature (20 °C -25 °C) before use. If crystals have formed in buffer

6 / 10

US and Canada: **Tel**: +1 800-810-0816

Web: http://www.acrobiosystems.com

E-mail: order@acrobiosystems.com




solution, warm to room temperature until the crystals have completely dissolved.

5. The kit should be stored at 2°C to 8°C.

TYPICAL DATA

The following data is for reference only. The sample concentration was calculated based on the results of the standard curve.

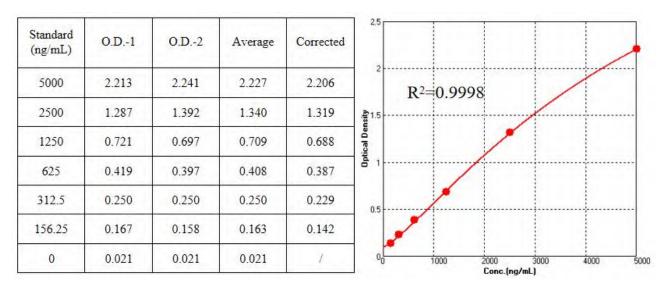

PRECISION

1. Intra-assay Precision

Three samples of known concentration were tested ten times on one plate to assess intra-assay precision.

2. Inter-assay Precision

Three samples of known concentration were tested in three separate assays to assess inter-assay precision.

	I	ntra-assay Precision	n	I	nter-assay Precision	n
Sample	1	2	3	1	2	3
n	10	10	10	3	3	3
Mean (ng/mL)	3896.658	933.040	400.426	4090.889	918.876	411.623
SD	104.589	49.678	20.014	168.258	45.840	20.479
CV (%)	2.7	5.3	5.0	4.1	5.0	5.0

Note: The example data is for reference only.

7 / 10

US and Canada:

Tel: +1 800-810-0816

Web: http://www.acrobiosystems.com

Asia and Pacific: Tel: +86 400-682-2521

E-mail: order@acrobiosystems.com





RECOVERY

Three samples with different concentrations were tested to calculate the recovery rate.

Sample(n=5)	Average Recovery %	Range %
High	91.8	81.8-113.6
Middle	91.5	85.5-100.5
Low	112.6	103.8-117.4

LINEARITY

To assess the linearity of the assay, samples spiked with high concentrations were serially diluted with calibrator diluent to produce samples with values within the dynamic range of the assay.

		Cell culture medium (DMEM)	Cell culture medium (1640)	
1.2	Average Recovery (%)	97.1	100.3	
1:2	Range (%)	94.3-99.9	98.1-102.2	
1:4	Average Recovery (%)	92.5	95.6	
1.4	Range (%)	91.3-94.6	92.7-97.6	
1:8	Average Recovery (%)	101.1	102.5	
1.0	Range (%)	100.5-101.5	100.1-103.9	
1:16	Average Recovery (%)	116.1	113.7	
1.10	Range (%)	113.0-118.9	111.8-116.2	

Note: The example data is for reference only.





PLATE LAYOUT

	1	2	3	4	5	6	7	8	9	10	11	12
Α	Std1	Std1		()						()	()	()
В	Std2	Std2		($$ $)$					(()
С	Std3	Std3	$\left(\begin{array}{c} \cdots \\ \end{array} \right)$	$(\underline{\cdot})$	\bigcirc		($$ $)$	($$ $)$	()	()	($$ $)$	()
D	Std4	Std4	(($)$	=		($)$	((₎	()	($ $	()
E	Std5	Std5	((\mathbb{Z}	\bigcirc	$\langle \cdots \rangle$	(()	$\langle z \rangle$	()
F	Std6	Std6	()	(\cdots)	$\stackrel{\dots}{\searrow}$	=	()	$\langle \cdots \rangle$		<u></u>	()	()
G	Blank	Blank	<u></u>		$\stackrel{\dots}{\nearrow}$	$\stackrel{\dots}{\searrow}$	$\langle \cdots \rangle$	()	;;; ;;;;;;;;;;;;;;;;;;;;;;;;;;;;;;;;;;	<u></u>	()	()
Н	Blank	Blank	()	()	()	()	()	()	("	()	()	()

Note: Blank is a Blank Dilution Buffer hole.

TROUBLESHOOTING GUIDE

Problem	Cause	Solution
Poor standard curve	* Inaccurate pipetting	* Check pipettes
Lauga CV	* Inaccurate pipetting	* Check pipettes
Large CV	* Air bubbles in wells	* Remove bubbles in wells
Walle day of	* Plate is insufficiently washed	* Review the manual for proper wash.
High background	* Contaminated wash buffer	* Make fresh wash buffer
Very low readings across the	* Incorrect wavelengths	* Check filters/reader
plate	* Insufficient development time	* Increase development time



RA195-EN.01

Samples are reading too high, but standard curve looks fine	* Samples contain cytokine levels above assay range	* Dilute samples and run again
		* Assay set-up should be continuous - have all standards and samples prepared appropriately before commencement
Drift	* Interrupted assay set-up	of theassay
	* Reagents not at room temperature	* Ensure that all reagents are at room temperature before
		pipetting into the wells unless otherwise instructed in the
		antibody inserts